📀 Trans Genic Inc.

Code No.KT042

For research use only

## Anti Rat MARCKS Monoclonal Antibody (Clone No. MAR11/2)

Myristoylated alanine-rich C kinase substrate (MARCKS), a specific substrate for protein kinase C, is abundant in the brain. MARCKS binds calmodulin under  $Ca^{2+}$  or actin under non- $Ca^{2+}$  conditions. Phosphorylation of MARCKS (P-MARCKS) can be used as an indicator of protein kinase C activation in intact cells. It has been shown that P-MARCKS has relatively weaker binding activity to calmodulin and actin compared with MARCKS.


Pathophysiologically, P-MARCS is also detected in the brain of Alzheimer's disease patients and in the smooth muscle cells of arteriosclerosis.

This antibody is specific to MARCKS and could be used for western blotting and immunohistchemistry.

| Package Size        | 50µg (200µL/vial)   |
|---------------------|---|
| Format              | Rat monoclonal antibody 0.25 mg/mL  |
| Buffer              | Block Ace as a stabilizer, containing 0.1% Proclin as a bacteriostat  |
| Storage             | Store below $-20^{\circ}$ C   |
|                     | Once thawed, store at $4^{\circ}$ C. Repeated freeze-thaw cycles should be avoided.   |
| Clone No.           | MAR11/2   |
| Subclass            | IgG2a   |
| Purification method | The lymphocytes from rat, immunized with C terminal peptides of rat MARCKS, were fused to mouse myeloma SP2 cells. The cell line (MAR11/2) with positive reaction was grown in serum free medium from which the |
|                     | antibody was purified by Protein G affinity chromatography.   |

Working dilution for immunohistochemistry: about 1~ 5µg/mL

The reactivity of anti MARCKS antibody using western blotting.

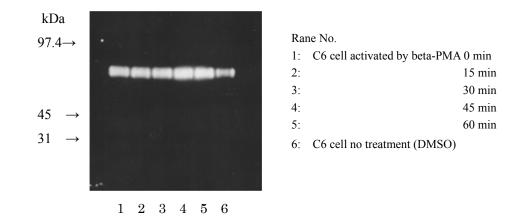

## Anti Rat MARCKS Monoclonal Antibody (Clone No. MAR11/2)

## [Reference]

 SD Rose et al, Differential expression of MARCKS and other calmodulin binding protein kinase C substrate in cultured neuroblastoma and glinoma cells.
J.Neurochem.63, 2314-2323,1994

Manufacturer



7-1-14 Minatojimaminami-machi, Chuo-ku, Kobe, Japan 650-0047 Telephone: +81-78-306-0295 URL:http://www.transgenic.co.jp techstaff@transgenic.co.jp